


# Form 2

Bring Dental 3D Printing to Your Desktop

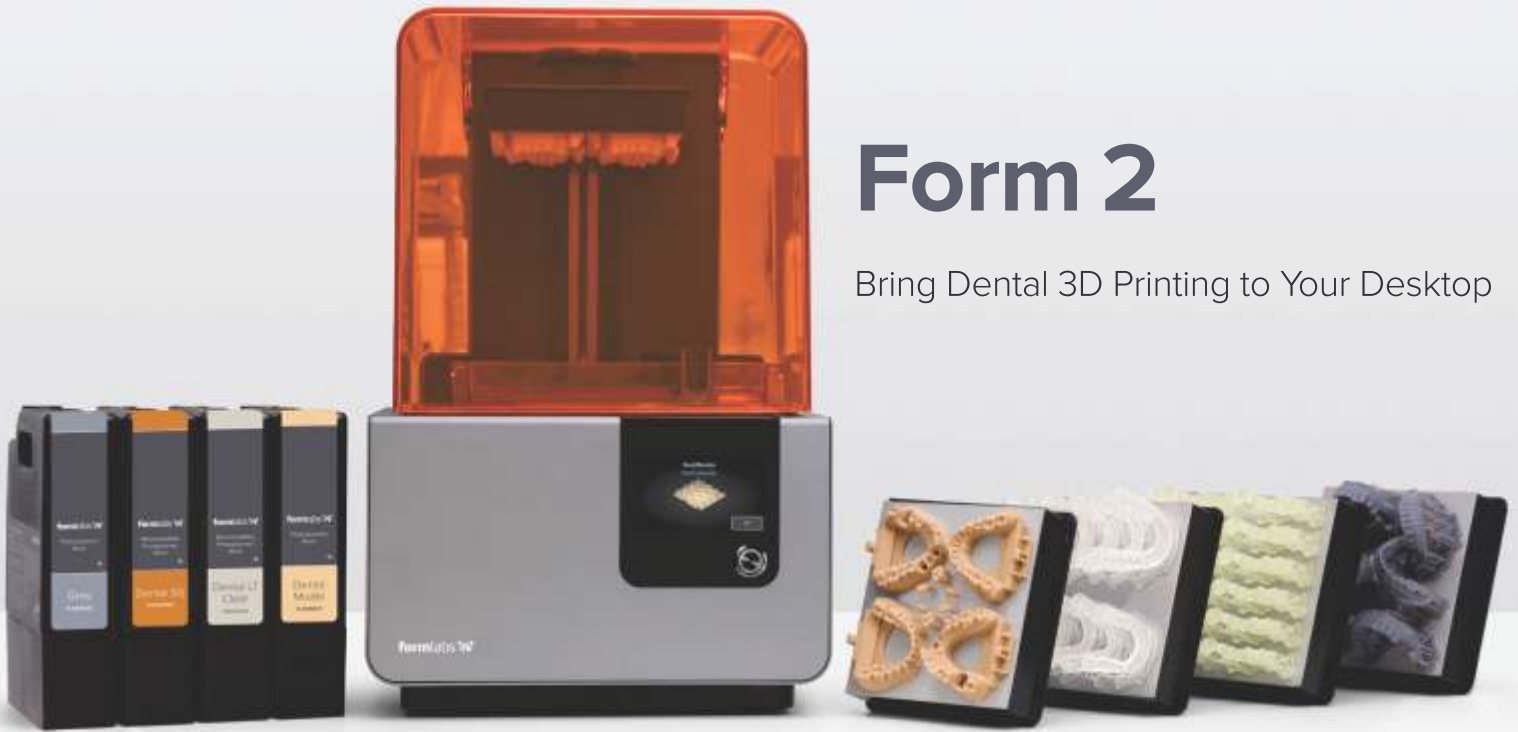

The Form 2 makes Dental 3D printing affordable for dental businesses of all sizes, enabling cost-effective, mass-customized digital production of a wide range of dental products.

## Fast Return on Investment

Affordable materials and low-cost machine pricing make 3D printing the cost-effective digital manufacturing choice for dental businesses of any size. Get a return on investment within months, rather than years.

Printing on the Form 2	Quad with Dies (Dental Model)	Full Arch with Dies (Dental Model)	Quarter Arch Surgical Guide	Full Arch Surgical Guide	Splint / Retainer (Dental LT Clear)	Full Arch Orthodontic Model (Grey)
PARTS PER BUILD	16 quads and dies	4 arches and dies	18 guides	12 guides	12 splints	9 arches
TIME AT 50 MICRON LAYER HEIGHT	2 quads / 3 hr 16 quads / 8.5 hrs	2 arches / 5.5 hr 4 arches / 8 hr	1 guide / 2.5 hr 18 guides / 6 hr	1 guide / 2.5 hr 8 guides / 7.5 hr	1 splint / 2.5 hr 12 splints / 9.5 hr	1 arch / 2.5 hr 9 arches / 7.5 hr
COST PER PART	₹ 110 / Quad	₹ 225 - ₹ 375/ Arch	₹ 145 - ₹ 225/ Guide	₹ 225 - ₹ 375/ Guide	₹ 295 - ₹ 450/ Split	₹ 110 - ₹ 265/ Split

# Dental Resins

With a library of dental materials, a single printer can be used for a wide range of applications. Automated resin cartridges and interchangeable resin tanks can be swapped in less than a minute, with no mess or wasted material.



### STANDARD GREY RESIN

Main Applications : Orthodontic, diagnostic, and educational models.



### DENTAL MODEL RESIN

Main Applications: Crown and bridge models with removable dies.



### DENTAL LT CLEAR RESIN

Main Applications: Splints, retainers, and other direct-printed orthodontic devices.




### DENTAL SG RESIN

Main Applications: Biocompatible Surgical guides and pilot drill guides.

## Printing Essentials

Streamline your workflows with the right accessories for consistent, successful printing

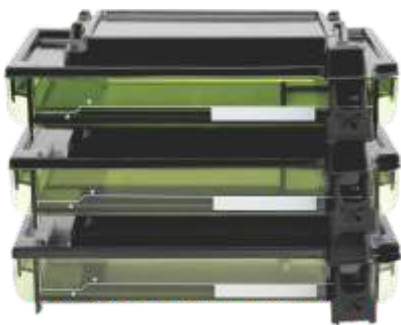

### Resin Tank LT

Long term tank for extended lifetime & high volume performance.



### Form Wash

Automated washing for consistently clean parts.



### Form Cure

Post-curing tuned and tested for each of our Dental Resins.

**UNICORN**  
**DenMart**   
Behind The Best Indian Dentists