


- ✦ Independent up & down movement & backrest movement of the chair which can be adjusted by the doctor as per convenience
- ✦ Double articulated headrest for handicapped & pediatric patients
- ✦ Removable dentist's & patient's feet system
- ✦ Automatic right arm rotatable at auto return to zero, for easy entry & exit of the patient
- ✦ 7 programs in the chair : 4 programmable working positions, other are spitting, last working position & auto return to zero
- ✦ With synchronised movement of back & seat (Trendelburg Position)
- ✦ Chair with anticrushing & movement locking system
- ✦ Soft comfort massage system for seat, back or both simultaneously
- ✦ Integrated bio-system for disinfection of handpiece tubings internally after every treatment
- ✦ Doctor's stool with adjustable backrest tilt & height adjustable, foot ring for attaining a correct working posture by the dentist especially during long dental treatment
- ✦ 3 axis movement of Cold & White 5 LED touchless sensor light (35,000 Lux) with variable intensity
- ✦ Removable & autoclavable handles for dental light head
- ✦ With original monitor arm

MULTIFUNCTIONAL FOOT CONTROL

- ✦ Emergency red switch to lock the chair in working position in case of any accidental obstruction
- ✦ One switch for 4 working position, one switch for last working / spitting position, one switch for auto return to zero, two independent switches for seat & backrest movement
- ✦ One switch for reversal of micromotor, one switch for light on / off
- ✦ Joy stick foot pedal for activation of air rotor handpiece with chip blower / water cut facility

ASSISTANT ARM

- ✦ Chromotherapy switch
- ✦ Water heating system in 3 way syringe
- ✦ Emergency stop on / off
- ✦ 2 points - one for low vacuum suction (saliva) and other for high vacuum suction cannulas only (blood + other debris)
- ✦ Water activation for spittoon & cup filler

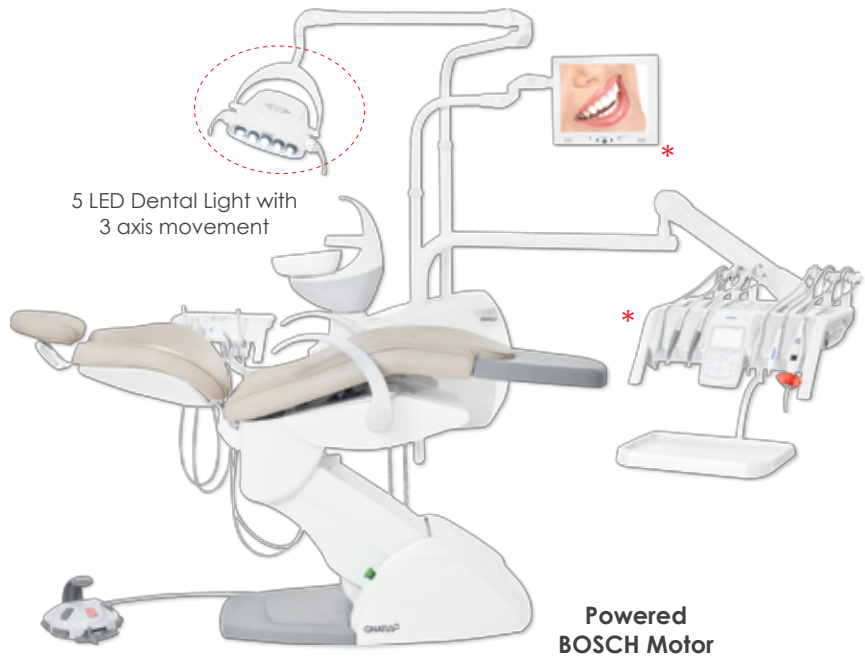

CE ISO 13485:2012 + AC:2012 / ISO 9001:2008

OVERHEAD DELIVERY UNIT

- ✦ With Pneumatic locking arm to prevent any unnecessary displacement of the Unit
- ✦ Feather touch control panel for easy touch play operations
- ✦ 3 points for air rotor with option for fibre optic handpiece
- ✦ 1 point for inbuilt airmotor/ micromotor
- ✦ An option for inbuilt scaler for oral prophylaxis
- ✦ One 3 way syringe
- ✦ Removable & autoclavable silicone pads for complete germ free operations
- ✦ A pressure gauze on the side of the delivery unit
- ✦ With Inbuilt LED X-Ray viewer

WATER UNIT

- ✦ 180° movable water unit
- ✦ Removable spittoon for easy cleaning
- ✦ Programmable water flow for convenience of operations by the operator
- ✦ Chromotherapy system for stimulating comfort & well being

ZERO BACK ACHE STOOL



- ✦ Available with all Gnatus Chairs. Helps prevent Spondylitis & low back pain. GNATUS has developed with considerable in R&D investment
- ✦ Up & down soft movements
- ✦ Adjustments of height & tilt are made through a side lever
- ✦ Adjustable from 50 to 70 cm of height with an adjustable ring
- ✦ The backrest with tilt & height adjustment which provides comfort to the lumbar area