

**DENTAL CHAIR**

- ✦ Independent up & down and backrest movement of the chair
- ✦ Luxury 6 LED sensor based light 10,000 - 25,000 Lux
- ✦ Safety switch for water-air - electricity with single on/off switch control
- ✦ Programable Dental Chair

**DELIVERY UNIT**

- ✦ With pneumatic locking arm to prevent any unnecessary displacement of the unit
- ✦ Feather touch control panel for easy touch play operation
- ✦ 3 Air rotor points and extra empty/ vacant point
- ✦ One 3- Way syringe.
- ✦ Detachable and easy clean handpiece holder block
- ✦ Control panel with 13 keys for chair control and other functions
- ✦ Inbuilt LED X-Ray viewer

**WATER UNIT**

- ✦ Rotary Glass cuspidor
- ✦ Detachable cuspidor, easy to clean
- ✦ Water heater for luke warm water

**MULTIFUNCTION FOOT CONTROL**

- ✦ Unique metal handle
- ✦ Chair position control
- ✦ Bowl -- control
- ✦ Cup filling control
- ✦ 2 Independent switches for air & water
- ✦ Stainless steel tubing for foot control operation



#P01  
Light Blue



#P02  
Apple Green



#P03  
Light Green



#P04  
Diamond Blue



#P05  
Yellow




#P06  
Orange



**FULL FUNCTION CONTROL PANEL WITH 13 BUTTONS**

1. Setting
2. Cuspidor flush
3. Cup filling
4. Water heating
5. Zero position
6. Emergency position control (-5 degree)
7. Spitting position
8. Operation light control
9. X-Ray viewer control

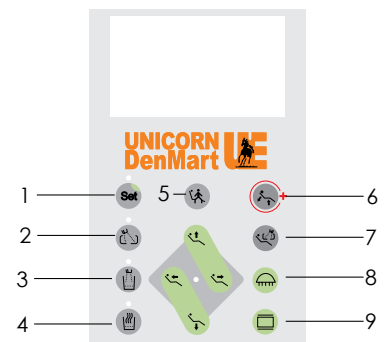

**ASSISTANT ARM**

- ✦ Low vacuum and high vacuum cannulas
- ✦ 3-Way syringe with option for luke warm water
- ✦ Chair control for assistant

**SPITTING POSITION SETTING**

On pressing the spitting position button, automatically the backrest moves upward, bowl rinsing start, light turns off, if pressed again, chair would automatically return to working status, backrest backward, bowl rinsing closes, lights turns on