

The DPEX III Gold Apex Locator is equipped with a clear bright LCD which gives you a crystal clear projection with different colours to clearly show the trajectory of the file. With the DPEX III Gold you get absolute precision leading to high RCT success rates. And what's more, with Unicorn Denmart after sales service, you have complete peace of mind.


- ✦ Based on advanced multiple frequency network impedance measurement technology and automatic calibrating ensures the measurements are accurate
- ✦ Autoclavable file clip, file hook, file probe. Avoiding cross infection effectively
- ✦ When doing the measurement, the changes of the sound can indicate the position of the file in the root-canal

PEAK PRECISION TO ENSURE PEAK RESULTS




1 Year Warranty


5 Years Warranty on PCB

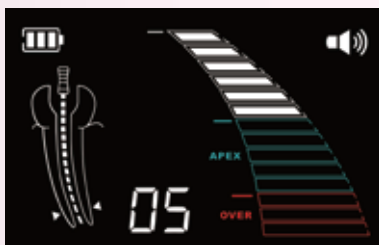

When the file reaches the front region of the apical foramen, the black bars displayed

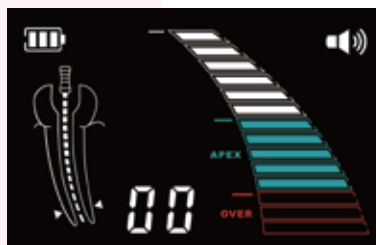

When the file reaches the position nearby the apical foramen, the green bars displayed

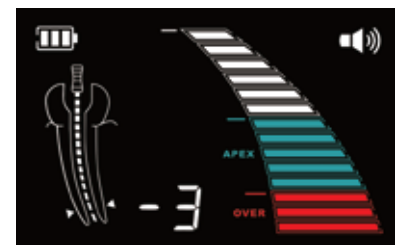

When the file exceed the apical foramen, the red bars display with a continuous beep sound

CE ISO 13485:2012 / AC:2012