

DENTAL CHAIR

- ✦ Independent up & down and backrest movement of the chair
- ✦ Luxury 6 LED sensor based light 10,000 - 25,000 Lux
- ✦ Safety switch for water-air - electricity with single on/off switch control
- ✦ Programmable Dental Chair

DELIVERY UNIT

- ✦ With pneumatic locking arm to prevent any unnecessary displacement of the unit
- ✦ Feather touch control panel for easy touch play operation
- ✦ 3 Air rotor points and extra empty/ vacant point
- ✦ One 3- Way syringe.
- ✦ Detachable and easy clean handpiece holder block
- ✦ Control panel with 13 keys for chair control and other functions
- ✦ Inbuilt LED X-Ray viewer
- ✦ Inbuilt woodpecker scaler

WATER UNIT

- ✦ Rotary Glass cuspidor
- ✦ Detachable cuspidor, easy to clean
- ✦ Water heater for lukewarm water


MULTIFUNCTION FOOT CONTROL

- ✦ Unique metal handle
- ✦ Chair position control
- ✦ Bowl - control
- ✦ Cup filling control
- ✦ 2 Independent switches for air & water



FULL FUNCTION CONTROL PANEL WITH 13 BUTTONS

1. Setting
2. Cuspidor flush
3. Cup filling
4. Water heating
5. Zero position
6. Emergency position control (-5 degree)
7. Spitting position
8. Operation light control
9. X-Ray viewer control

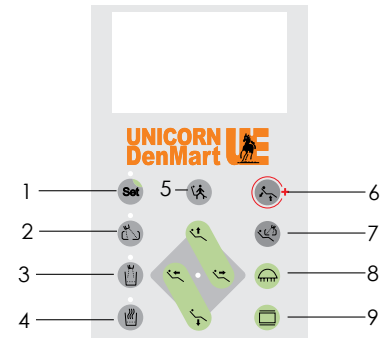

ASSISTANT ARM

- ✦ Low vacuum and high vacuum cannulas
- ✦ 3-Way syringe with option for lukewarm water
- ✦ Chair control for assistant

SPITTING POSITION SETTING

On pressing the spitting position button, automatically the backrest moves upward, bowl rinsing start, light turns off, if pressed again, chair would automatically return to working status, backrest backward, bowl rinsing closes, lights turns on