

## Surgic AP

Non Optic

### Main Features

#### Control Unit

- ✦ Max. Output/Max. Torque : 210 W / 50 Ncm
- ✦ Power Supply : AC 230 V 50/60 Hz
- ✦ Max. Pump Output : 75 mL/min
- ✦ Programs : 8 Programs
- ✦ Dimensions : W 275 x D 207 x H 101 (mm)

#### Micromotor (SG50MS)


- ✦ Motor Speed : 200-40,000 min<sup>-1</sup>
- ✦ Non-optic
- ✦ Brushless

### 8 Individual Program Options

NSK Surgic AP can memorise 8 complete program options. Each program includes speed, torque, coolant flow, and rotation direction.

### Large User-Friendly LCD

The Large LCD display offers immediate recognition of the selected operational program data.

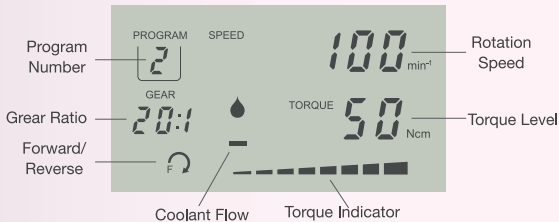

### User Friendly Foot Control

The Surgic AP foot control (FC-55) is simple to operate and conforms to IPX8 - the standard for medical foot control devices and certification for operating theatre use. (forward direction only)



## Provides Solid Performance for Surgical & Implant Procedures



With S-Max SG20 Handpiece

### Powerful External Pump

Surgic AP guarantees you optimum irrigation capacity and simple mounting with an external powerful water pump. The Surgic AP pump is easily detachable and can be repeatedly autoclaved.

